

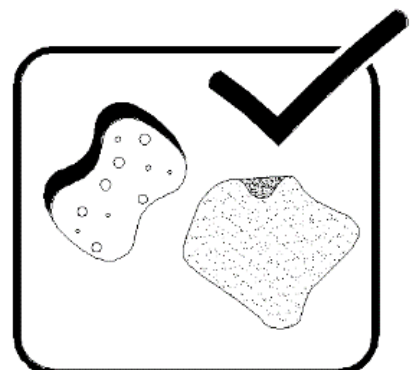
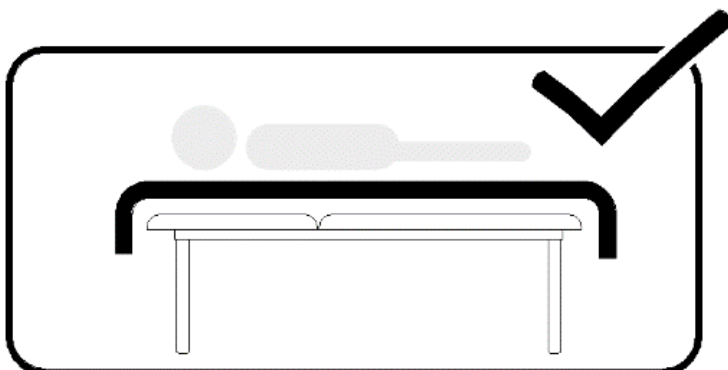


UPHOLSTERY

USEFUL INFORMATION AND RECOMMENDATIONS

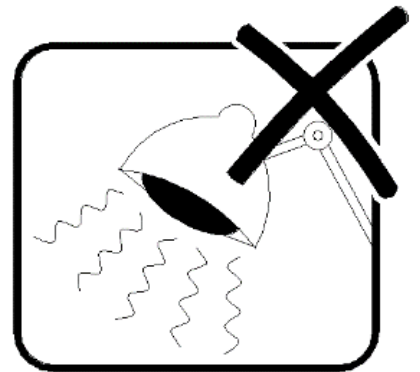

IMPORTANT

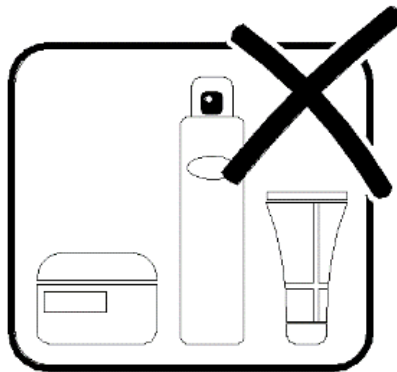

UPHOLSTERY

USEFUL INFORMATION AND RECOMMENDATIONS

Ecopostural's superior vinyl used in the manufacturing of this product has been supplied to us by Benecke-Kaliko (Grupo Continental) and meets the D60-070 standard governing household, food and pharmaceutical products.

Our vinyl is inherently flame retardant M2 class according to EC compliance. This means that only the area in contact with an ignition source or an incandescent material will consume but in a gradual and predominantly flameless way.

All our coated fabrics have the Sanitized® additive offering an antimicrobial and antibacterial protection. The vinyl technical features have been approved to offer excellent results.

In order to ensure the durability and keep a long lasting performance of our vinyl please follow our **RECOMMENDATIONS**:

Easy maintenance

- Do the dusting regularly with a white damp cloth.
- In the event of stubborn stains, clean using soapy water (a solution of water with a few drops of soap) and a smooth sponge. Then dry with a cloth.
- The use of different cleaning products other than the above-mentioned may damage the outer varnish coating aimed at protecting the vinyl against external agents.
- Avoid excess water when rinsing the vinyl as it may get into the seams.

Stain removing

- In the event of a stain, clean immediately using soapy water to avoid the stain persistence.
- Please refrain from using alcohol, all solvents and stain removers as they contain harmful components which will damage the vinyl.
- If you finally decide to use these products under your responsibility, it is highly recommended to previously test their results in a hidden area of the upholstery.

Prevention is better than cure

- The upholstery surfaces more specifically exposed to a direct contact with patients (face, chest, hands areas) are more likely to deteriorate due to the effect of perspiration and other types of exudation. Massage creams and oils may migrate into the coated fabric affecting its durability.
- The use of inappropriate products during maintenance may damage the upholstery.
- Certain clothing (i.e. denim) dyes may run onto lighter materials and migrate into the coated fabric mass.
- The use of infrared lamps and other heat sources may dry out and crack the vinyl.
- Be careful when carrying out specific techniques requiring the use of extremely heat sources (i.e. moxibustion)

The most significant traces of a damaged vinyl are the bright appearance and the rough touch. These signs will be preceding vinyl crack.

For these reasons we strongly recommend you to **PROTECT** your Ecopostural table or any other upholstered product with a thick cover sheet or even better with a protective towel in order to retain the original properties of the vinyl for a long time.

GUARANTEE

Benecke-Kaliko (Grupo Continental) cannot be held responsible in the event of an inappropriate use or maintenance which results in damage to the protective varnish on the coating leaving it brittle. ECOPOSTURAL is merely informing you and reminds you that our guarantee will not be applicable in those cases.